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Read this manual carefully before use, and retain it for future reference.

Safety Instructions


When using the toothbrush for the first time, it is common for the gums to experience itching or slight bleeding, which can be attributed to the fact that your gums have not yet adapted to the efficient cleaning provided by the toothbrush. It is recommended to start with the gentle cleaning mode. The conditions will typically be resolved within a period of one to two weeks. In the event that they persist for more than two weeks, please stop using the toothbrush immediately and seek medical assistance.

WARNING:

- The plug must be removed from the socket-outlet before cleaning or maintaining the appliance.
- This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.
- Children should be supervised to ensure that they do not play with the appliance.
- Do not expose the appliance or battery to excessive temperatures.
- Be aware of the risk of terminals of the battery-operated appliance or battery being short circuited by metal objects.
- Only charge the appliance at a temperature between 5°C and 35°C.
- This appliance can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved.
- Children shall not play with the appliance.
- Cleaning and user maintenance shall not be made by children without supervision.
- This appliance contains batteries that are non-replaceable.
- For the purposes of recharging the battery, only use the detachable supply unit provided with this appliance.

Note: Illustrations of product, accessories, and user interface in the user manual are for reference purposes only. Actual product and functions may vary due to product enhancements.

Product Overview

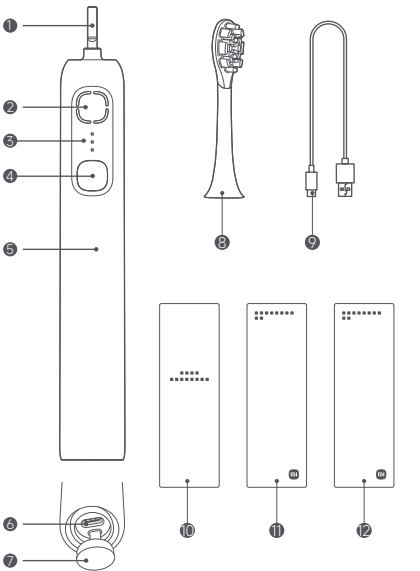- ① Shaft

③ Mode/Battery Level Indicator
 - Gentle Cleaning Mode
 - ⋮ Standard Cleaning Mode
 - ⋮ Deep Oscillation Mode

⑧ Toothbrush Head

⑩ Quick Start Guide

⑫ Warranty Notice
- ② Brushing Area Indicator

④ Power Button

⑤ Toothbrush Handle

⑥ Charging Port


⑦ Charging Port Cover

⑨ Type-C Charging Cable

⑪ User Manual

Installing the Toothbrush Head

Keep the bristle direction of the toothbrush head in the same direction as the front of the toothbrush handle, align downwards to install the toothbrush head on the handle, and make sure there is a slight gap between the toothbrush head and the handle (The gap is designed to prevent the toothbrush head from colliding with the handle during high-frequency vibration, so as to generate more power from the motor).

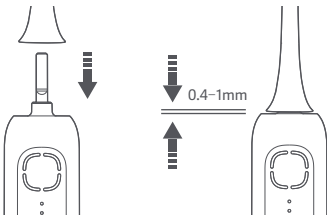How to Use

Button operations

Turning on: Press the button to turn on the toothbrush.

Switch modes/pause: Within 10 seconds after turning on the toothbrush, press the button to cycle through gentle cleaning, standard cleaning, deep oscillation modes and pause. After the toothbrush has been started for 10 seconds, press the button to pause it.

Turning off: After pausing or brushing, the toothbrush will turn off after 12 seconds of inactivity.

Pick-up-to-wake:

When the toothbrush is turned off, lifting it manually will wake it up, and the mode indicator will display the mode used before turning off or low battery status.

Note: If the wake-up feature is triggered multiple times in a short period, it will be temporarily disabled. The feature can be reactivated by either charging or turning the toothbrush on.

Introduction to modes

- Gentle cleaning mode: Suitable for people with sensitive teeth
- Standard cleaning mode: Suitable for daily oral hygiene
- Deep oscillation mode: Suitable for deep cleaning of stubborn tooth stains

Notes:



When the toothbrush is turned on, there is a 3-second gradual power increase function, and the vibration will gradually increase within 3 seconds.

Each brushing mode lasts for 2 minutes, after which the toothbrush will automatically pause. It is recommended to start with the gentle cleaning mode for new users.

Healthy brushing

For better brushing effect, it is recommended to use the Bass brushing method, gently moving the toothbrush with the bristles positioned at a 45° angle towards the gum line.

CAUTION: The toothbrush pauses briefly every 30 seconds to prompt the user to change the brushing area.

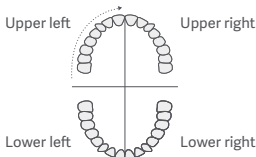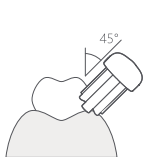30-second Area Pause

Missing area reminder

The brushing area indicator will provide feedback on the current brushing situation when brushing is paused or completed. The brushed area will illuminate in white, while the missing area will illuminate in orange.

CAUTION:

1. When brushing the missing area, after completing each area, there will be a brief pause of vibration to indicate the need to change the brushing area.
2. The duration of supplemental brushing varies depending on the areas missed during brushing. Once all the missed areas are brushed, the toothbrush will automatically pause and display the brushing results.
3. To achieve relatively accurate brushing results, it is recommended to maintain a stable head position and keep your gaze fixed straight ahead during the brushing process.

Orange: Missing area

White: Brushed area

Care & Maintenance

Cleaning the toothbrush head

1. After brushing teeth, rinse the toothbrush head with running water.
 2. After rinsing, shake off water and place the toothbrush head in a dry and well-ventilated area to allow for air drying.
- * Dentists recommend replacing the toothbrush head every three months.

Charging

Charging method

Open the charging port cover at the bottom of the handle and use the charging cable supplied with the toothbrush for charging.

Display of battery level indicator

Toothbrush Status	Battery Level	Battery Level Indicator Status
Operating	20%–100%	Displaying the mode, instead of battery level
	5%–20%	Blinking orange three times and then going off
The toothbrush cannot be turned on	0%–5%	Blinking orange; the toothbrush vibrates three times
Charging	100%	Keeping white
	0%–100%	Breathing white

Caution:

- It takes about six hours to get the toothbrush fully charged once the battery is depleted.
- Please use a 5 V $\overline{\sim}$ 1 A power adapter (sold separately) when charging.

Moisture detection of charging port

When connecting the charging cable to the toothbrush, if the battery level indicator continues to blink orange, it indicates the presence of liquid at the charging port. Charging at this time poses a risk of short-circuiting and please ensure that the charging port is dry before charging.

Unplug the charging cable

If there is any liquid present at the charging port, proceed with the following cleaning steps:

1. Unplug the charging cable.
2. Keep the charging port downwards and shake gently to drain excess liquid.
3. Place the toothbrush and the charging cable in a dry and well-ventilated place for air drying. Once the toothbrush is completely dry, try to recharge it. If the battery level indicator continues blinking orange, it is recommended to extend the drying period until normal charging resumes.

CAUTION: Please do not use paper towels, cotton swabs, etc. to dry the charging port to avoid foreign objects from entering.

Environmental Protection


When the toothbrush reaches the end of its service life, please remove, dispose of and recycle the built-in battery according to local laws and regulations of the country or region where the product is used. If the battery cannot be removed, dispose of it in strict accordance with the regulations on classified domestic waste disposal in your area.

⚠️ WARNINGS:


- Do not expose the battery to fire to prevent the battery from exploding or releasing toxic and hazardous substances.
- Do not short-circuit the battery to prevent it from combusting and causing a fire.
- Do not disassemble the toothbrush unless you want to dispose of it; otherwise, the toothbrush will be damaged.
- Before disposing of the toothbrush, remove the battery first.
- The battery must be removed from the appliance before it is scrapped.
- The appliance must be disconnected from the supply mains when removing the battery.
- The battery is to be disposed of safely.

Remove rechargeable battery


1. Insert a flat screwdriver into the gap between the handle's bottom cover and the toothbrush body, and pry upward until the bottom cover opens.

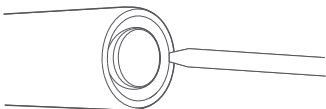2. Use a flat screwdriver to pry the buckles inwards on both sides of the inner compartment, press the front power shaft against a hard object, and push the inner compartment out of the toothbrush.

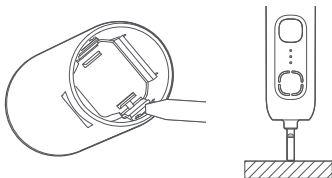3. Use a diagonal cutter to cut off the battery connection wires. Remove the battery and deposit it at an official recycling point.

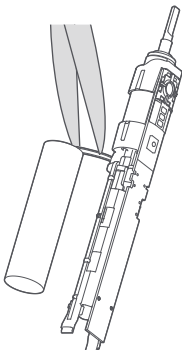Troubleshooting

Issues	Possible causes	Solutions
The toothbrush does not turn on.	The toothbrush is charging.	Unplug the charging cable from the toothbrush.
	The battery level is too low.	Charge the battery in a timely manner.
The toothbrush head cannot be installed.	There are foreign objects inside the toothbrush head.	Remove foreign objects.
	There are foreign objects on the shaft.	Remove foreign objects.
	Non-genuine toothbrush head is used.	Replace with a genuine toothbrush head.
The toothbrush head vibration is weak.	The toothbrush head is not installed in place.	Reinstall the toothbrush head.
	There is something stuck between the toothbrush head and the handle.	Detach the toothbrush head and remove the foreign objects.
The battery level indicator blinks orange six times and the toothbrush cannot be turned on.	The operating temperature is too high or too low.	Please use at the ambient temperature from 5°C to 35°C.

Specifications

Name: Electric Toothbrush	Model: MES609
Rated Power: 2.5 W	Rated Input: 5V --- 1A
Charging Time: Approx. 6 hours	IP Rating: IPX8

WEEE Disposal and Recycling Information

All products bearing this symbol are waste electrical and electronic equipment (WEEE as in directive 2012/19/EU) which should not be mixed with unsorted household waste. Instead, you should protect human health and the environment by handing over your waste equipment to a designated collection point for the recycling of waste electrical and electronic equipment, appointed by the government or local authorities. Correct disposal and recycling will help prevent potential negative consequences to the environment and human health. Please contact the installer or local authorities for more information about the location as well as terms and conditions of such collection points.

Manufacturer: Xiaomi Communications Co., Ltd.
Address: #019, 9th Floor, Building 6, 33 Xi'erqi Middle Road, Haidian District, Beijing, 100085, China
For further information, please go to www.mi.com
User Manual Version: V1.0
Made in China

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製造元：Xiaomi Communications Co., Ltd.

所在地：#019, 9th Floor, Building 6, 33 Xi'erqi Middle Road, Haidian District, Beijing, 100085, China

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WARRANTY NOTICE

As a Xiaomi consumer, you benefit under certain conditions from additional warranties. Xiaomi offers specific consumer warranty benefits which are in addition to, and not instead of, any legal warranties provided by your national consumer law. The duration and conditions related to the legal warranties are provided by respective local laws. For more information about the consumer warranty benefits, please refer to Xiaomi's official website

<https://www.mi.com/global/support/warranty>. Except as prohibited by laws or otherwise promised by Xiaomi, the after-sales services shall be limited to the country or region of the original purchase. Under the consumer warranty, to the fullest extent permitted by law, Xiaomi will, at its discretion, repair, replace or refund your product. Normal wear and tear, force majeure, abuse or damage caused by the user's negligence or fault are not warranted. The contact person for the after-sale service may be any person in Xiaomi's authorized service network, Xiaomi's authorized distributors or the final vendor who sold the products to you. If in doubt please contact the relevant person as Xiaomi may identify.

The present warranties do not apply in Hong Kong and Taiwan.

Products which were not duly imported and/or were not duly manufactured by Xiaomi and/or were not duly acquired from Xiaomi or a Xiaomi's official seller are not covered by the present warranties. As per applicable law you may benefit from warranties from the non-official retailer who sold the product. Therefore, Xiaomi invites you to contact the retailer from whom you purchased the product.

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