Have a glow up with Spark

User Manual

Soocas Spark Sonic Electric Toothbrush


Table of Contents

User Manual (EN) -----	07
User Manual (JP) -----	14
User Manual (KR) -----	21
User Manual (RU) -----	28

Thank you for choosing Soocas oral care series: Soocas Spark Sonic Electric Toothbrush

Please read this manual carefully and the correct method of usage for Soocas Spark Sonic Electric Toothbrush to fully understand its performance to better clean and care for your oral health.

Product Name: Soocas Spark Sonic Electric Toothbrush

Product Model: MT1

Rated Input: 5V 0.5A

Charging Time: 2h

Charging Interface: Type-C

Store Temperature: 0-45°C/32-113°F, Humidity

Made in China

Manufacturer: Shenzhen Soocas Technology Co., Ltd.

Address: Room 1101, 1102, 1103, 1104, 401 & 402, Building 2, Chongwen Area, Nanshan iPark, No. 3370 Liuxian Avenue, Nanshan District, Shenzhen, China

Table of Contents

• Precautions	01
• Quick Start	03
• Product Description.....	03
• User Manual	04
• Proper Brushing Techniques	05
• Properly Using the Soocas Spark Sonic Electric Toothbrush	06
• Environmental Notice	07

Precautions

In order to use this product in the correct manner to avoid any personal injury or property loss, please read all instructions and warnings carefully before using the product.

Warning

1. This appliance can be used by children aged 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning the use of the appliance in a safe way and understand the hazards.
2. Do not charge in a wet environment.
3. Do not use this product while charging.
4. This appliance is not allowed for children to play.
5. Users who have orthodontic correction, tooth repair, or implant teeth should use this product under the guidance of dentists.
6. This product complies with the safety standards of electro-magnetic devices. For any problems related to cardiac pacemaker or other implant devices, please consult the doctor or the manufacturer of the implant device before us
7. This product is only designed to clean the teeth, gums, and tongue. Do not use it for any other purpose. If you feel any discomfort or pain, please stop using this product and contact the doctor.

8. If the power cord is damaged, it must be replaced with a special cord, or a special component purchased from its manufacturer or maintenance department.
9. Before charging, please ensure that the charging port of this product is dry.
10. Do not use this product when the charging port cover is not properly placed.
11. This appliance must only be supplied at SELV.

Caution


1. Do not place the charging cable equipped with this product in a wet environment or water.
2. This product contains a built-in rechargeable battery. Do not put it into the fire for heating, or charge, use, or place it at high temperature.
3. Do not modify, disassemble, or repair this product by yourself.
4. When charging, please use a 5 V charger (purchased separately).
5. Do not use paint thinner, benzene, alcohol, or other chemicals to clean this product.
6. Do not charge this product for a long time without care. The maximum charging time should not exceed 12 hours.
7. If this product is not used for a long time, please fully charge it at least once every six months.

Important

Before using the product, please check the package for any obvious damages. Please do not use the product if the package is damaged.

Quick Start

1. Product Description

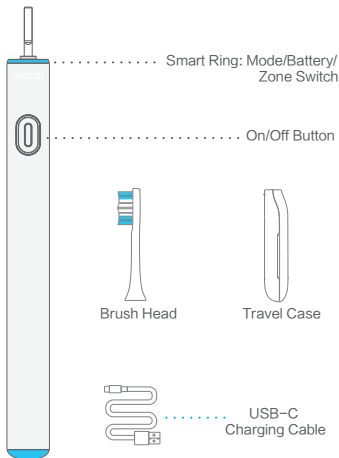2. User Manual

- Initiation of the Spark

- Short press to turn on. After 10 seconds, short press to turn off.
- Press consecutively within 10 seconds for the followings: turn on, switch modes, turn off.
- Long press for longer than 1 second: Activate pause mode. Within 30 seconds of the pause mode, the timer will be paused.

- Gentle Mode

Gentle vibrations to clean teeth. Recommended for first-time users and those with sensitive teeth.

- Cleaning Mode

High-frequency vibration for effective cleaning. Ideal for whitening of teeth.

- Charging & Battery Status Indicator

- Flashing Red Light: Low battery, please charge at your earliest convenience.
- Flashing Blue Light: Charging indicator, the product is currently under charging.
- Constant Blue Light: The product has been fully charged.


- Zone Switch Reminder

Every 30 seconds, the Smart Ring flashes with a slight vibration pause as a reminder to switch brushing zones.


3. Proper Brushing Techniques

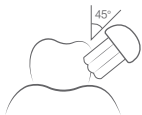1. Guide the brush head slowly from tooth to tooth with a 40-degree angle between the bristles and the gum lines.

2. Press the On/Off button for the preferred brushing mode.

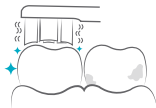3. The toothbrush will vibrate to indicate the user has to brush a new area of the mouth.

4. Keep the brush head to the teeth closely without bending the bristles.

Properly Using the Soocas Spark Sonic Electric Toothbrush

Battery Status and Charging

1. Soocas Spark Sonic Electric Toothbrush has a built-in rechargeable battery with a working duration of up to 40 days after it is fully charged.
2. Soocas Spark Sonic Electric Toothbrush has been fully charged and discharged for aging before leaving the factory. However, there may be a deviation with the declared standby time upon the initial use due to natural discharge in the production period. However, you can still immediately enjoy your extreme cleaning experience with the product after unpacking.
3. Soocas Spark Sonic Electric Toothbrush is equipped with a USB-C high-efficiency charging port, which features high compatibility and fast charging speed. You can charge the toothbrush using any charging adapter in compliance with national standards, and please keep the charging port dry and tightly sealed. However, using a poor-quality charger may bring unpredictable security risks.
4. According to the test, it takes about 2 hours to fully charge the battery. The specific charging time may be subject to the remaining power and the charging environment.

Environmental Notice

When the product has reached the end of its service life, you may need to remove the built-in battery for recycling before discarding the product according to the national or local laws and regulations.

Warning:

Please do not throw the battery into a fire to prevent the battery from exploding or releasing toxic and hazardous substances. Please do not make the battery short-circuited to prevent the battery from burning and causing a fire.

Unless you would like to discard the product, do not disassemble the product yourself, as premature disassembly will damage the product.

Take basic safety preventive measures before operations. The power supply must be disconnected before removing the battery.

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The first part of the paper discusses the importance of understanding the underlying mechanisms of the observed phenomena. This involves a thorough review of the existing literature and a clear identification of the research gaps. The second part presents the methodology used in the study, which includes a combination of qualitative and quantitative approaches. The data collection process is described in detail, along with the statistical analysis performed. The results of the study are then presented, showing the key findings and their implications. Finally, the paper concludes with a discussion of the limitations of the study and suggestions for future research.

The study aims to explore the relationship between the independent variable and the dependent variable, taking into account the moderating effect of the third variable. The research hypothesis is that there is a positive correlation between the independent variable and the dependent variable, which is moderated by the third variable. The study is designed to test this hypothesis using a sample of participants from various backgrounds.

The data was collected through a series of surveys and interviews, ensuring a comprehensive understanding of the participants' experiences. The analysis revealed that the independent variable has a significant positive impact on the dependent variable, and this relationship is further strengthened by the presence of the third variable. These findings have important implications for the field of study and provide valuable insights for practitioners.

In conclusion, the study successfully tested the research hypothesis and provided evidence for the proposed model. The findings suggest that the independent variable is a key factor in explaining the variation in the dependent variable, and the third variable plays a crucial role in moderating this relationship. Future research should continue to explore these relationships in different contexts and with larger samples to further validate the findings.

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